



Agriculture and  
Agri-Food Canada

Agriculture et  
Agroalimentaire Canada



# Algae Identification

## LAB GUIDE

Agriculture and Agri-Food Canada  
Agri-Environment Services Branch

---

© Her Majesty the Queen in  
Right of Canada, 2011  
Cat. No. A125-8/1-2011E-PDF  
ISBN 978-1-100-18307-7  
AAFC No. 11432E

---

Aussi offert en français sous le titre :  
*Identification des algues :*  
*Guide de Laboratoire*

---

For additional copies of this publication or to  
request an alternate format, please contact:  
Publications Services

**Agriculture and Agri-Food Canada**

Tower 4, Floor 0, Room 160  
1341 Baseline Road  
Ottawa (Ontario) K1A 0C5 Canada  
Telephone: (613) 773-1444  
Fax: (613) 773-1498  
E-mail: [publications@agr.gc.ca](mailto:publications@agr.gc.ca)

---

**Disclaimer**

Responsibility for interpretation of the content  
contained in this publication rests with the user.  
Information in this publication is provided solely  
for the user's information and, while thought to be  
accurate, is provided strictly "as is" and without  
warranty of any kind, either expressed or implied.  
Agriculture and Agri-Food Canada will not be liable  
to you for any damages, direct or indirect, or lost  
profits arising out of your use of this publication.

---

# Algae Identification **Lab Guide**

ACCOMPANYING MANUAL TO THE  
ALGAE IDENTIFICATION FIELD GUIDE

Agriculture and Agri-Food Canada  
Agri-Environment Services Branch

Nancy Serediak  
Streamline Consulting  
Edmonton, AB Canada



Mai-Linh Huynh  
Agriculture and Agri-Food Canada  
Regina, SK Canada

## Acknowledgments

The authors sincerely thank the reviewers and supporters of this project: AAFC/AESB Water Quality Division; Glen Brandt, AAFC-AESB; Garth Mottershead, AAFC-AESB; Wayne Wark, AAFC Communications; Mark Graham, University of Alberta; Ron Zurawell, Government of Alberta; and Bob Klemmer, Government of Saskatchewan. Their comments and recommendations have proven invaluable for ensuring the relevance and scientific accuracy of the reference manual and field guide. The authors are also grateful to the following for allowing use of their photos: Charles Delwiche, David Patterson (micro\*scope website), Morgan Vis, David Krogmann and Mark Schneegurt (Cyanosite website), Mark Graham, Steve Murrell, Ron Zurawell, David John, Brian Whitton, Peter York, Jane Jamieson, and Nick Stewart.

The authors would also like to acknowledge Steve Murrell, AAFC/AESB Water Quality Division, for coordinating and leading the quality assurance component of this project.

The authors would also like to sincerely thank Michael Parry, AAFC-AESB GIS Unit in Calgary Alberta, for his hard work and collaboration in developing the GIS-based program tool for this project.

# Table of Contents

Understanding the Algal Groups	5-7
Taxonomic List	8-9
Blue-Greens	10-18
Green Algae	19-32
Diatoms	33-35
Flagellates	36-38
Other	39
Glossary	41
Photo credits	42-43
References	44
Recommended Reading	45
On-line Resources and Identification Keys	46

# Understanding the algal groups

Understanding the common algal groups is the first step to identifying algae. Algae are classified into four major groups (Palmer 1962): Blue-greens, Greens, Diatoms and Flagellates.

By knowing the major group of your algae, verification and identification to the genus and species level with the use of a microscope and classification key becomes easier. Many traditional classification keys can be complex and daunting if you do not have an idea where to start. It is our hope that this lab reference manual will provide a good starting point.

Below are characteristics that each Algal Group will typically exhibit, however be advised that exceptions for some species do exist.

## Blue Greens (Known as cyanobacteria)

- ▶ Includes species that fall under the Division Cyanophyta, Kingdom Eubacteria.
- ▶ They are referred to as algae but are actually photosynthetic bacteria that appear like true algae.
- ▶ They are mostly single-celled organisms or aggregations of single-celled organisms (planktonic). Planktonic algae generally do not adhere together in a mass, in that you could not easily grab a handful or mass of these organisms. Some species, such as Nostoc, are gelatinous masses that can be picked up.
- ▶ Anabaena, Aphanizomenon, Oscillatoria and Myrocystis often dominate plankton communities and are notorious for producing potentially toxic blooms in fresh waters (note: only certain species produce toxins). These toxins can cause death and illness in livestock; and gastroenteritis, liver damage and skin and eye irritations in humans.
- ▶ High numbers of blue-greens often indicate high nutrient levels in the water, although Oscillatoria can form dense accumulations in less fertile waters.
- ▶ \*Cell colour is uniform throughout (i.e. no membrane-bound coloured structures called chloroplasts) and may appear blue-green to violet, but sometimes red or green.
- ▶ \*They have very simple internal structure, no definite nucleus or intracellular membrane-bound organelles.

## Greens

- ▶ Includes species from the Division Chlorophyta (Green Algae).
- ▶ They are morphologically diverse (micro and macroscopically). Microscopically, they can be uni- or multi-cellular, colonial, have branched or unbranched filaments and vary in size.
- ▶ Macroscopically, they can be planktonic, filamentous or found attached to objects (epiphytic).

- ▶ Filamentous forms of algae grow as long multi-cellular threadlike filaments or masses (i.e. you can pick it up with your fingers). Early in the season these algae may develop clumps scattered on the water surface but as the season progresses they may cover the entire surface of a water body.
- ▶ \*Some have a whip-like tails called a flagella. Swimming cells have 2,4, or 8 flagellae of equal length at the anterior end.
- ▶ \*All species have internal membrane-bound structures (chloroplasts) that are often bright green. Other structures within the cell may also be seen at times, especially starch granules which stain dark with a preservative know as Lugols iodine solution.

\* Characteristics only visible through a microscope.

## Diatoms

- ▶ Includes species from the Class Bacillariophyta, of the Division Chrysophyta.
- ▶ They are often too small to see with the naked eye.
- ▶ Many are single celled, and less often colonial or filamentous. They may be planktonic, although many grow on other submerged objects.
- ▶ \*They are easily identified by the cell wall (composed of silica) and shape, which are often complex and ornately marked.
- ▶ \*They have shades of green, yellow, or brownish coloured chloroplasts.

## Flagellates

- ▶ Includes species from various Divisions (including Chrysophyta, Euglenophyta, and Pyrrophyta) that are motile, i.e. they have whip-like tails called flagellae.
- ▶ \*They are often microscopic uni-cellular (solitary) planktonic algae.
- ▶ \*They have shades of golden brown to yellow or green in the chloroplasts, or can even be colourless. Colourlessness may indicate that the species is totally dependant on het erotrophic nutrition. Some have a hard external covering (lorica) which makes them appear brown or black.

\* Characteristics only visible through a microscope.

*Images in this manual are microscopic images, unless otherwise indicated.  
For some images, magnification or scale was not available.*

---

# Taxonomic List\*

of algal genera referenced in this manual

## Phylum (Division) Chlorophyta (Domain Eukaryota; Kingdom Protista)

Sub-phylum Chlorophyceae	
Order Volvocales	
Family Chlamydomonadaceae	
Chlamydomonas	pg 19
Family Volvocaceae	
Volvox	pg 20
Order Tetrasporales	
Family Tetrasporaceae	
Tetraspora	pg 21
Order Chlorococcales	
Family Oocystaceae	
Ankistrodesmus	pg 22
Chlorella	pg 23
Family Scenedesmaceae	
Scenedesmus	pg 24
Family Hydrodictyaceae	
Hydrodictyon	pg 25
Pediastrum	pg 26
Order Chaetophorales	
Family Chaetophoraceae	
Stigeoclonium	pg 27
Order Oedogoniales	
Family Oedogoniaceae	
Oedogonium	pg 28
Order Siphonodadales (Cladophorales)	
Family Cladophoraceae	
Cladophora	pg 29
Order Zygnematales	
Family Zygnemataceae	
Spirogyra	pg 30
Family Desmidiaceae	
Closterium	pg 31
Sub-phylum Charophyceae	
Order Charales	
Family Characeae	
Chara	pg 32

### Phylum (Division) Euglenophyta (Domain Eukaryota; Kingdom Protista)

Order Euglenales	
<i>Euglena</i>	pg 36

\*Adapted arrangement from G.W. Prescott (1978).

### Phylum Chrysophyta (Domain Eukaryota; Kingdom Protista)

Sub-phylum Chrysophyceae	
Order Ochromomadales	
Family Dinobryaceae	
<i>Dinobryon</i>	pg 38
Family Synuraceae	
<i>Synura</i>	pg 37
Sub-phylum Bacillariophyceae	
Order Pennales	
Family Fragilariaceae	
<i>Asterionella</i>	pg 33
<i>Synedra</i>	pg 34
Order Centrales	
Family Coscinodiscaceae	
<i>Cyclotella</i>	pg 33

### Phylum Cyanophyta (Domain Eubacteria, Kingdom Monera)

Order Chroococcales	
Family Chroococcaceae	
<i>Anacyctis/Microcystis</i>	pg 10
<i>Chroococcus</i>	pg 11
<i>Gloeocapsa</i>	pg 12
Oder Oscillatoriaceae	
Family Oscillatoriaceae	
<i>Oscillatoria</i>	pg 13
<i>Spirulina</i>	pg 14
Order Nostocales	
Family Nostocaceae	
<i>Anabaena</i>	pg 15
<i>Aphanizomenon</i>	pg 18
<i>Nostoc</i>	pg 16
Family Rivulariaceae	
<i>Gloeotrichia</i>	pg 17

# BG Blue Greens

## ► MICROCYSTIS

*Microcystis* and *Anabaena* (pg15) blooms appear very similar. Heavy blooms of these species create the memorable pea-soup conditions associated with cyanobacteria. Colours range from green to blue-green. Although large colonies are readily visible to the naked eye, individual cells are not visible to the naked eye and identification to genus or species requires a compound microscope. The most common species is *M. aeruginosa*; *M. flos-aquae* and *M. wesenbergi* can also co-occur with *M. aeruginosa*.

*Microcystis* is found throughout the water column because like other blue-greens, they possess gas vesicles (small air-filled structures) that allow the colonies to regulate buoyancy. This adaptation allows *Microcystis* to remain at optimal light levels to attain maximum photosynthetic rate.


Blooms occur after periods of hot, calm weather and can develop rapidly over a few hours. This is because during periods of wind induced mixing of the water column, *Microcystis* have to produce more gas vesicles to counter the downward drag of the current. But when the wind ceases, the cells are unable to reduce the number of gas vesicles and as a result, they will migrate en masse to the surface over a few hours.

Most large blooms occur from mid to late-summer, and can through to fall.

*Microcystis* species are capable of producing microcystin, a hepatotoxin (liver toxin) that is lethal to wildlife, pets, livestock and humans that consume untreated or improperly treated water. Heavy blooms may cause skin irritation for recreational water users (dermatitis), as compounds on the surface of *Microcystis* colonies are skin irritants.

**Under the microscope:** *Microcystis* colonies are better defined and more densely clustered relative to an *Anabaena* colony. Cells are spherical shaped (~4-5µm in diameter), enclosed in a mucilaginous sheath and usually joined with other individuals within the sheath to form small detached colonies; one can find hundreds of cells in one sheath. Each cell has gas vesicles that often appear blackish blue. Cells are non-filamentous, planktonic and not attached to objects.


N.B. *Microcystis* is often used synonymously with *Anacystis* in some literature. A new naming approach, introduced by Francis Drouet (1968), has resulted in the disuse of the name *Anacystis*.

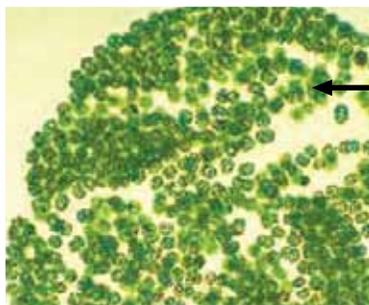

Dense colony enclosed in mucilaginous sheath

Spherical shaped cells

4-5µm diameter

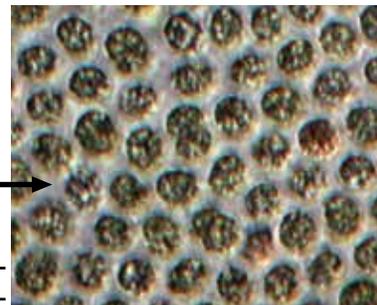

D. Krogmann and M. Schneegurt  
(Cyanosite 2006)


D. Patterson (micro\*scope 2006)

---

► **CHROOCOCCUS**

---

Similar to *Gleocapsa* (pg12), *Chroococcus* can be found attached to rocks or soil in or near water but have a colourless mucilaginous sheath. Under the microscope, cells are non-filamentous, planktonic and solitary or grouped in small families of 2, 4 or 8.

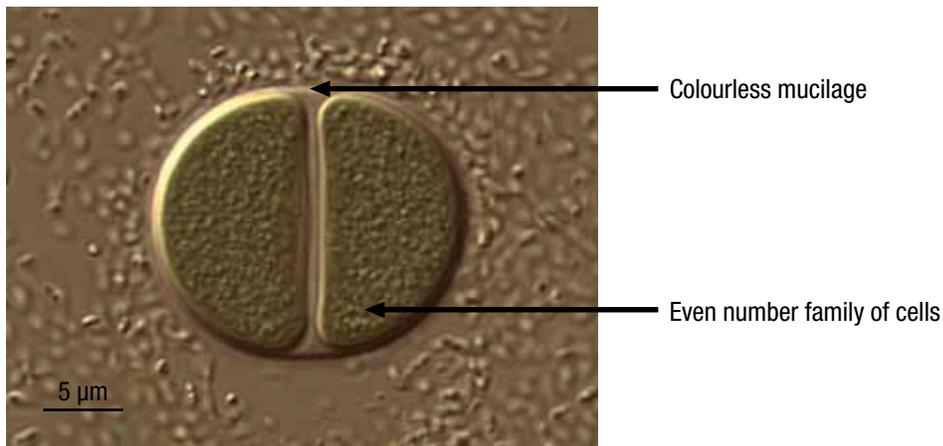


D. Patterson (micro\*scope 2006)

# BG Blue Greens

## ► GLOEOCAPSA

This subaerial genus often forms red, brown or orange gelatinous globs on substrates in or near water. The jelly-like colonies can be quite large and occur commonly on damp soil and rocks, but more typically on moist soil.

**Under the microscope:** When colonies are composed of only a few cells, *Gloeocapsa* can be confused microscopically with *Chroococcus* (pg11). These spherical non-flagellated, non-filamentous cells can be solitary or grouped in 2, 4, or 8 individuals that remain aggregated in gelatinous masses. Each individual cell is enclosed in a thick, homogenous or lamellate, mucilaginous sheath. The sheath may be yellow, brown or red, and even blue to violet.

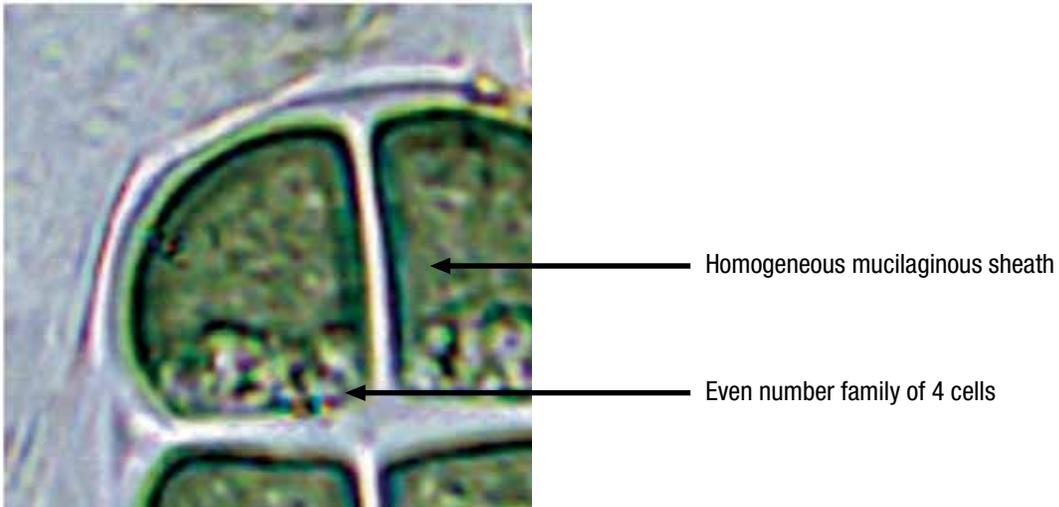

M. Vis (Ohio University 2002)

---

## ► OSCILLATORIA

---


This filamentous blue-green can produce microcystin (hepatotoxin) and anatoxin-a (neurotoxin). *Oscillatoria* do not bloom in highly fertile waters, but blooms can occur in less productive (oligotrophic) lakes and reservoirs.

*Oscillatoria* species often range in colour such as green to blue, blue-green, purple and red. Individuals can occur singly, in aggregates with other algal communities, free-floating or attached to an object. *Oscillatoria* colonies occurring on rocks are one of the reasons fishermen wear felt-soled boots.

*Oscillatoria* blooms tend to occur at depths **in the water column** other than at the water's surface. Thus, water may appear clear on the surface but there may be a bloom forming three metres down. This can be important for people who have buried water intake lines in a surface water source. Blooms can be mat-forming as well. If it is left in a shallow dish in the lab, it will creep up the sides or spread itself over the bottom.

*Oscillatoria rubescens* is known to produce blooms under ice, turning the water brilliant red to dark purple. These blooms usually occur during late winter but may persist several weeks following ice melt in spring and contribute a red colour to the open water.

**Under the microscope:** *Oscillatoria* trichomes (thread of cells) often appear as dense, opaque ribbons which oscillate under the microscope. Hence, their name comes from their characteristic wavy 'oscillating' motions. The cells are short cylinders, usually wider than long. Trichomes are straight or irregularly twisted, with a swollen apical cell.

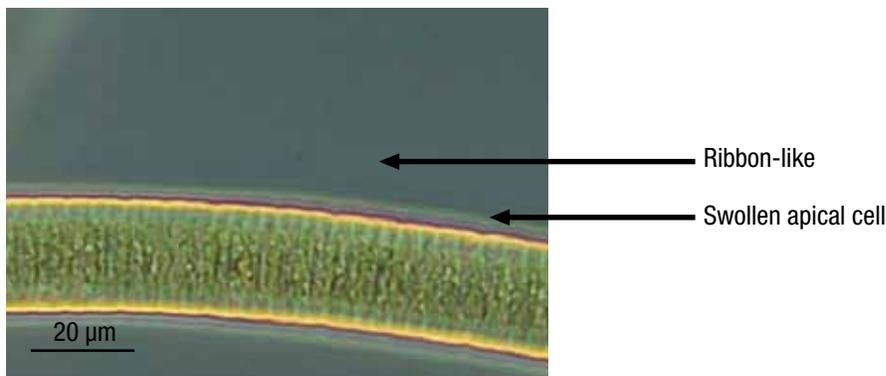

B. Anderson and D. Patterson (microscope 2006)


# BG Blue Greens

## ► SPIRULINA

This genus does not form the intense blooms often associated with cyanobacteria nor does it produce toxins, but is of particular interest to agricultural producers because of its market value. It is most commonly marketed by health food producers as a general curative.

In the field, *Spirulina* colonies are commonly found as somewhat amorphous, light green masses either by themselves or associated with *Oscillatoria* (pg13). They are evident in standing water bodies from late spring through to fall.

**Under the microscope:** *Spirulina* appears as a coiled thread of cells (trichomes) and are unicellular. Individual trichomes are essentially comprised of single cells that spiral down its entire length. Species are differentiated in part based on the tightness of the spiral. Like many blue-greens, *Spirulina* species extrude mucilage which flows along the length of the thread. This mucilage creates movement which can be viewed microscopically.

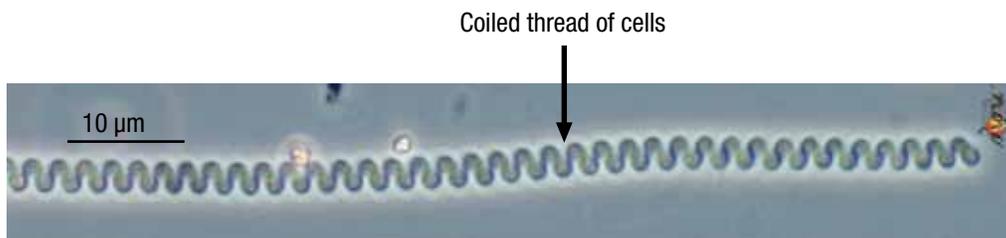

D. Patterson, L.A. Zettler and V. Edgecomb (micro\*scope 2006)

---

## ► ANABAENA

---


This genus is filamentous but does not form large colonies, and thus colonies are still distinguished in the field as planktonic with the 'open finger test'. Some *Anabaena* species are capable of producing anatoxin, or microcystin, or both. Although less common than *Microcystis*, periodic blooms of *Anabaena* can cause deaths of pets, livestock and wildlife. Heavy blooms may cause skin irritation for recreational water users.

*Anabaena* species are prevalent in phosphorus rich water and can cause pea-soup blooms in standing water after periods of warm, calm weather from late spring through summer. *Anabaena* can be especially prevalent if water is nutrient rich. There are often noxious taste and odour concerns associated with *Anabaena* blooms, especially once the bloom begins to decay. You may see tufts or flocs with aquatic plants or films on sediment.

**Under the microscope:** *Anabaena* is filamentous or threadlike, and characterized by the trichomes (thread of cells without its sheath) that are often coiled/spiralled. *Anabaena* threads often clump into loosely associated colonies held together by a soft mucilage. Trichomes are not parallel to each other, are multi-cellular, beadlike and not bundle forming. Each cell is spherical, beadlike, and non-flagellate.



M. Graham and S. Murrell (400X magnification)

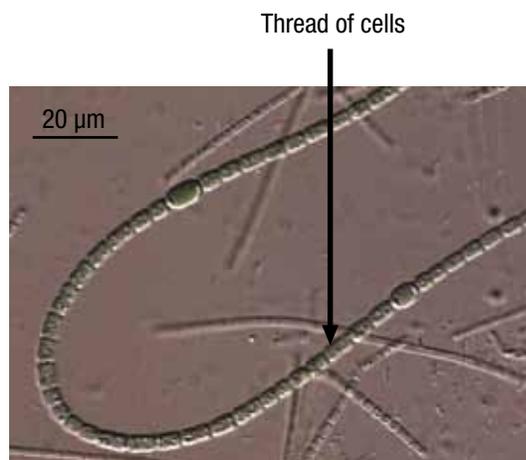

Australian Biological Resources Study (microscope 2006)


# BG Blue Greens

## ► NOSTOC

*Nostoc* is probably one of the most readily identified cyanobacteria in the field. The colonies are found in a characteristic jelly-like ball which is actually a thick mucilage surrounding heavily coiled, bead-like threads. The rubbery colonies are often found floating in standing water bodies from mid-summer through to fall. Care should be taken with young children swimming in surface waters in summer. *Nostoc* colonies have sometimes been called 'freshwater grapes', and children have been known to eat them. Since the colonies are capable of producing cyanotoxins, this is not a recommended practice.

Colonies rarely occur in great numbers. They can be macroscopic in size, growing up to a few centimetres in diameter (between a dime to quarter in size). Colour of the mucilage ranges from yellow-clear to somewhat blueish to olive green.

**Under the microscope:** A colony contains a tangled mass of trichomes (thread of cells) that are enclosed in a thick mucilaginous matrix. Cells of the trichome are bead-like or barrel-shaped similar to *Anabaena* (pg15), but the trichomes are randomly tangled in a common jelly-ball sheath.

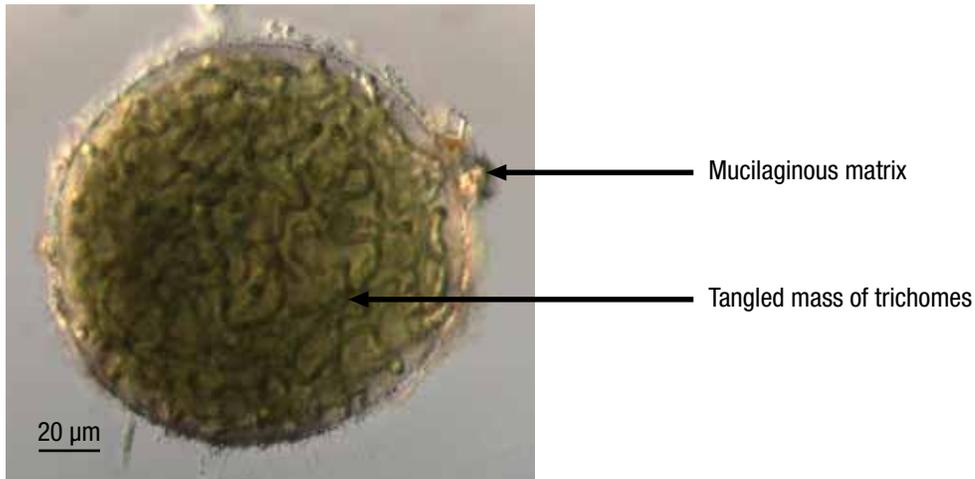

M. Bahr and D. Patterson (Micro\*scope 2006)

---


## ► GLOEOTRICHIA

---

*Gloeotrichia* is generally recognized as non-toxic, but there is some question as to whether they are capable of producing toxins. Heavy *Gloeotrichia* blooms can however cause serious skin irritations for recreational swimmers. When mature, colonies are visible to the naked eye upon close inspection and look like numerous, tiny fish eyes, tapioca beads or pom poms. A heavy bloom will make the water appear buff or tan in colour. Some species of this genus can be green or blue. *Gloeotrichia* occurs primarily in lakes, ponds and other standing water, and blooms can occur from mid-summer to early fall.

Some *Gloeotrichia* species are found attached to submerged aquatic plants. Colonies are of a similar shape, but the mucilage tends to be much firmer in these epiphytic species. Colour ranges from green to black in these forms, and colony size is from 1 to 2 millimeters in diameter. Planktonic species are more soft. Common species found is *G.echinulata* and can occur in abundance in the plankton of hard water lakes.

**Under the microscope:** Individual threads of cells (trichomes) are enclosed in a soft mucilage resulting in small, globular colonies. Trichomes are straight or irregularly twisted, and tapered from base to apex.

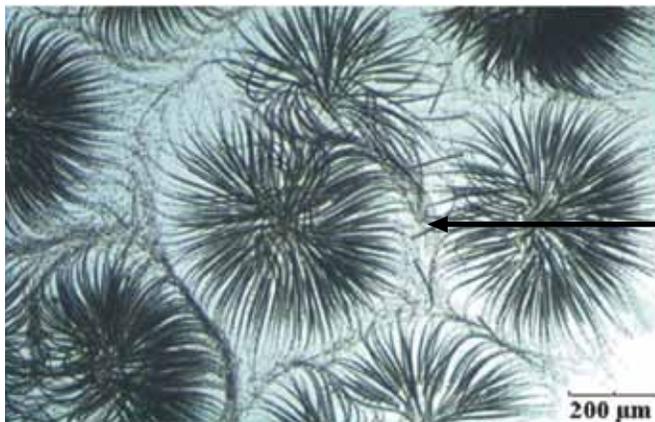

P. York (York et al. 2002)

## BG Blue Greens


### ▶ APHANIZOMENON

*Aphanizomenon* is one of the most easily recognized genus of the blue-greens. Individuals form threads that tend to parallel one another and join into colonies to form a distinctive shape that resembles tiny grass or green fingernail clippings. In natural settings, the colonies are sickle-shaped clusters roughly 0.5 to 1.0 centimeters long. If cultured in a lab, colonies will not readily form this shape. The most common species is *A. flos-aquae*.

Heavy *Aphanizomenon* blooms also create pea-soup conditions in the water column and form after extended periods of warm, calm weather from late spring through summer. When in bloom, they tend to join into loosely formed clumps that superficially resemble some green algal colonies. However, scooped handfuls will easily break apart and pass through the open fingers of a cupped hand.

Some species are capable of releasing a neurotoxin called saxitoxin, similar to the paralytic shellfish poison that occur in marine environments. In North American freshwater environments, *Aphanizomenon* appears to be primarily non-toxic though it can cause skin irritation to recreational water users. Decaying blooms also can release high concentrations of ammonia, potentially leading to summer fish kills.

**Under the microscope:** This genus is filamentous or threadlike, where the thread of cells (trichomes) are straight and parallel, forming a free-floating flake-like bundle. Each trichome contains 1-3 heterocysts. Heterocysts are large cells that look different in appearance- thick walled, weakly pigmented or clear compared to the other cells. Each trichome has tapered ends which appear like chopped grass. A distinctive characteristic of this genus is the akinete cell, which is cylindrical, darker in colour compared to the heterocysts and other cells, and typically found adjacent to the heterocysts. Akinete cells can help identify this genus to the species level

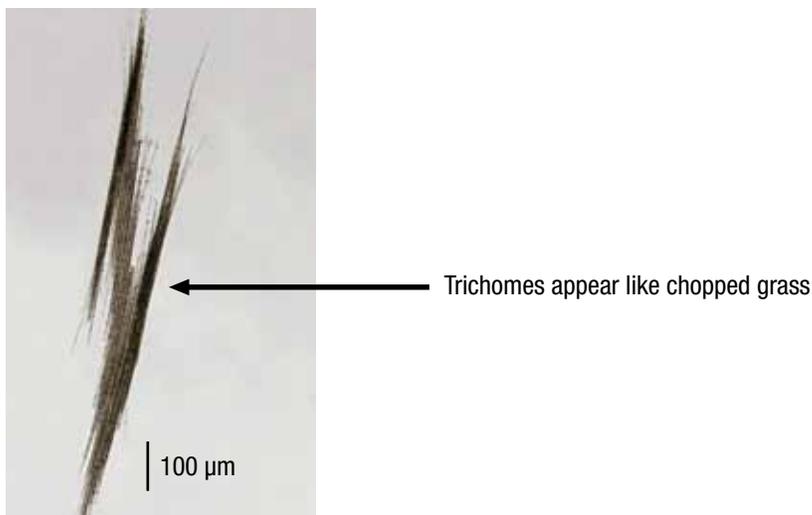

W. Bourland (micro\*scope 2006)


---

## ► CHLAMYDOMONAS

---

Species of this genus are found in a wide variety of habitats, even in snow fields in the high alpine. They are common bloomers that exhibit a fishy smell in high numbers.

**Under the microscope:** They are solitary, small, elliptical-oval cells. They are also biflagellate (i.e. they have two whip-like tails). A tiny red **eyespot** at the base of the flagella is often evident. There is also a single large cup-shaped **chloroplast** that occupies most of the volume of the cell. They are often easily confused with other biflagellated algae or with the flagellated reproductive stages of other genera.

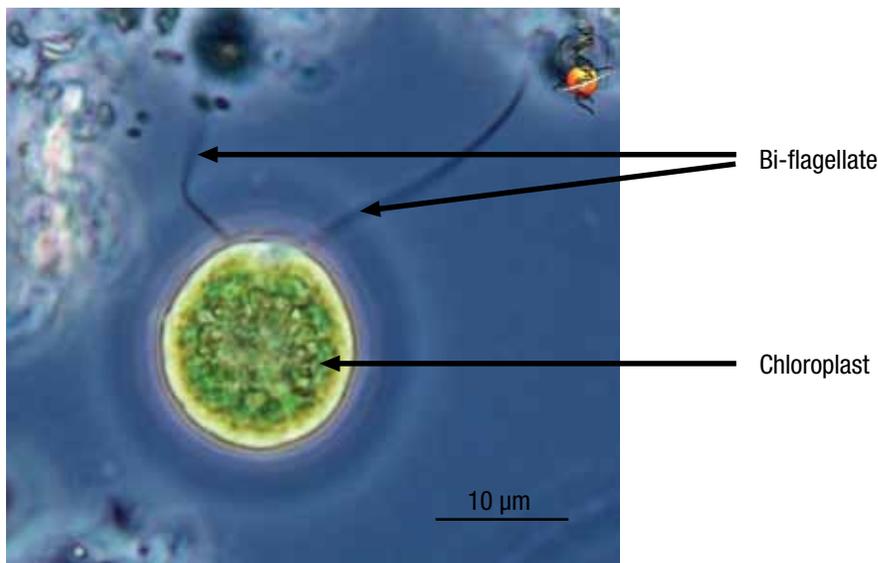



D. Patterson and A. Laderman (micro\*scope 2006)

# GR Green Algae

## ► VOLVOX

*Volvox* is often found within the general plankton of standing water systems. Blooms are infrequent but will occur over short periods during warm months, especially in waters contaminated with nitrogenous wastes or organically enriched water. Individuals of this genus join to form complex, spherical colonies that at their largest are visible to the naked eye (about 1 millimeter in diameter) and are usually yellowish-brown to light green.

**Under the microscope:** It is a large hollow sphere of a colony of algal cells, usually 500 or more, surrounded by a gelatinous matrix. Individuals comprising the sphere are spherical or ellipsoid, and biflagellate, i.e. they have two whip-like tails of equal length. The synchronized beating of all flagellae propels the colony through the water column. The cells are in a single layer on the periphery of the colony. Most cells are vegetative. Each cell also contains green cup-shaped internal structures called chloroplasts and large internal structures that store starch, which would stain dark with iodine.

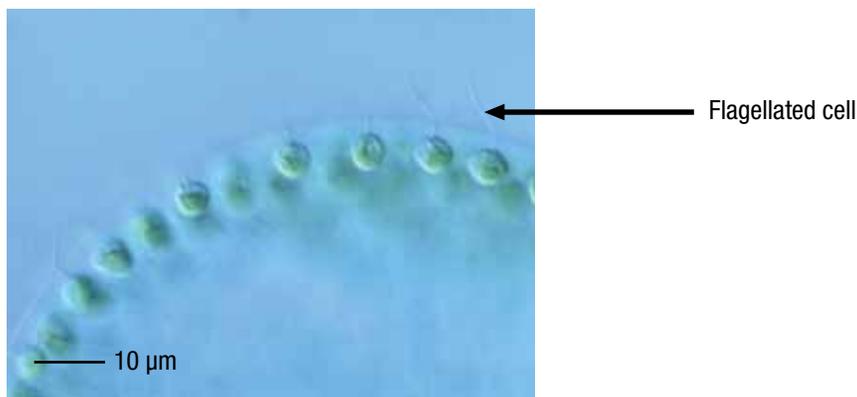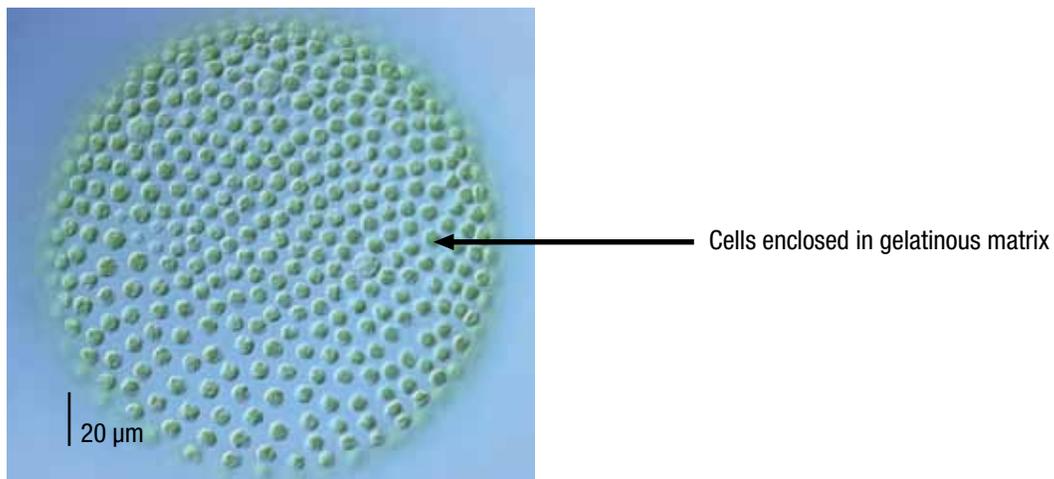


D. Patterson and M. Farmer (micro\*scope 2006)

---

► **TETRASPORA**

---

*Tetraspora* species have colonies of minute, rounded cells in a macroscopic soft, fragile, mucilaginous common tube or sac. They are one of the first to develop in early spring. They are non-motile, non-filamentous and epiphytic (i.e. found attached to aquatic plants). They are typically associated with submerged substrata or can be found floating as colonial clumps and massive sheet-like skins in ditches, ponds, bogs and lakes.

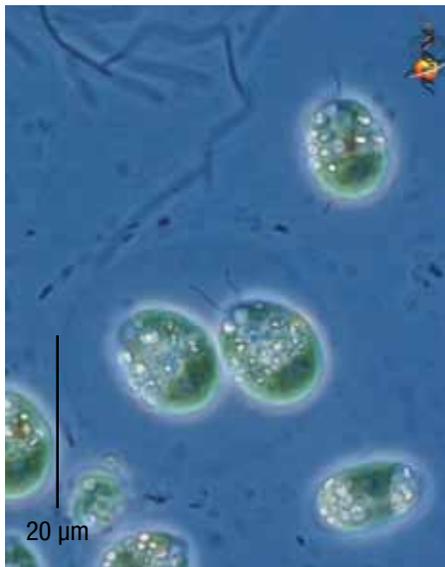


D. Patterson and M. Farmer (micro\*scope 2006)

## GR Green Algae

### ► ANKISTRODESMUS

These species are common within the open water phytoplankton. They will occasionally form pale green blooms **within the water column** during hotter months, but can bloom with regularity in personal aquaria or artificial pools. They are actually indicators of cleaner water since they tend to disappear from the algal community in very polluted systems.

**Under the microscope:** *Ankistrodesmus* exists as single cells which may form loose associations with each other (bundles, entangled with others), but not true colonies. Cells often elongate to cylindrical, straight, curved or spirally twisted shapes.

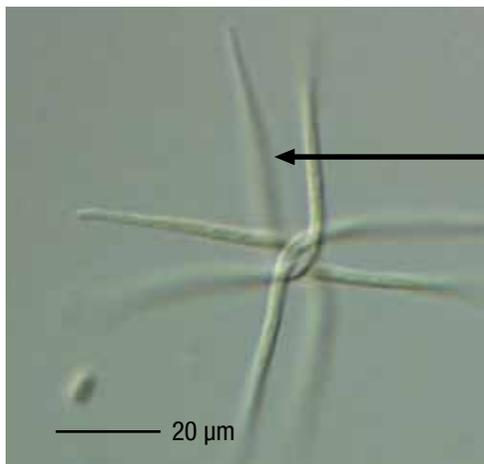

← Spirally twisted cell, entangled with another cell

M. Bahr and D. Patterson (micro\*scope 2006)

---

## ► CHLORELLA

---

*Chlorella* are green single-celled algae widely distributed in various fresh and marine water, soil and subaerial habitats. They are often consumed by freshwater invertebrate animals such as sponges. It has a grass-like smell because of the high amounts of chlorophyll contained within it, the highest concentration of any plant in the world. Certain species of *Chlorella* are harvested as a health food for humans and lower grades of *Chlorella* are fed to livestock. They are known to contain more than 20 vitamins and minerals including B complex, beta-carotene, vitamins C and E, iron, calcium, up to 70% protein, and 19 of the 22 amino acids.

**Under the microscope:** *Chlorella* is a small spherical green cell, often solitary and is non-flagellated. Its chloroplasts are thin and cup-shaped.




W. Bourland (micro\*scope 2006)

# GR Green Algae

## ► SCENEDESMUS

Colonies are too small to be visible to the naked eye. Blooms can occur in natural systems during warmer months, colouring the water a lovely shade of green. It is frequently abundant in nutrient-rich waters. This genus will typically develop along with others belonging to the Order *Chlorococcales* such as *Ankistrodesmus*, *Chlorella*, *Hydrodictyon* and *Pediastrum*. *Scenedesmus* is also a common uninvited guest of many freshwater aquaria. *Scenedesmus* is widely distributed in standing water systems. It occurs as a typical member of the plankton of open and near-shore water. The most common species is *S. quadricauda*.

**Under the microscope:** Colonies in this genus occur in a wide variety of forms but are comprised of 4-12 cells adjoined along the entire length of 1 cell wall. They form a coenobium, a colony that always has the same number of cells once it is formed. Cell arrangement and shape is used to distinguish one species from another, except in laboratory cultures where development is often atypical. Colonies contain chloroplasts, and are non-motile and non-filamentous. Cell walls are smooth, or with 1 or 2 curved spines or teeth. Note that not all species of this genus have the typical spines on the cell walls.

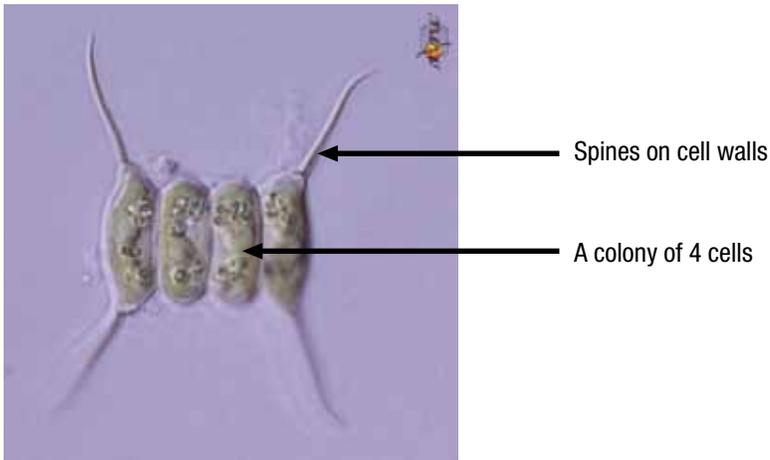

L.A. Zettler and D. Patterson (micro\*scope 2006)

---


## ► HYDRODICTYON

---

*Hydrodictyon* is a filamentous alga which forms net-like shaped colonies. The common name for this genus is 'water net'. This is owing to its habit of growing so profusely that it can become a very bothersome, weed-like intruder. It is common in summer in standing or very slow-moving water, including irrigation ditches. It does not produce toxins but it will bind boat props and clog intake lines. Scooping up a handful of the alga will reveal the net-like appearance of the colony. This genus is found growing in hard water lakes and ponds. Dense mats are unsightly in appearance and produce unpleasant odours particularly as the organism begins to decompose. This genus is also very susceptible to copper treatment.

**Under the microscope:** Cells can be seen grouped into cyclical or sheet-like colonies to form definite patterns, and a network of cells. Reproduction in each cell produces not just another cell, but a replication of the entire original structure. In this manner, *Hydrodictyon* grows very rapidly and mat sizes quickly become macroscopic. They are non-motile, non-filamentous and colonial.

A network of cells

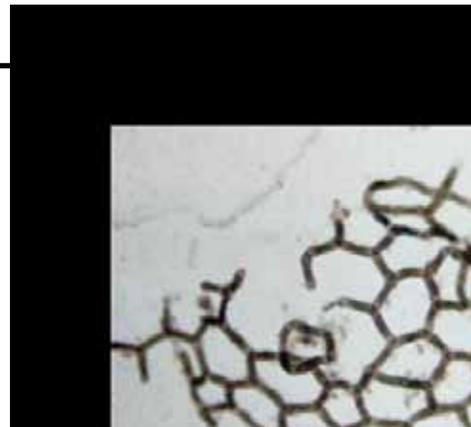


M. Vis (Ohio University 2002)

## GR Green Algae

### ► PEDIASTRUM

On occasion, *Pediastrum* species will form deep green blooms during hot summer months. Blooms are usually **relegated to the water column** and do not result in the surface scums associated with cyanobacteria (do not produce film on substrate). Colonies commonly occur in ponds, small lakes and ditches intermingled with other members of the plankton community. The cell walls of *Pediastrum* are so extraordinarily rigid and resistant to decay that specimens have been found in fossil records. Blooms are often common in moderately nutrient-enriched water bodies.

**Under the microscope:** *Pediastrum* is colonial with a very distinctive, unmistakable star shape. Colonies are flat, circular and have a multi-cellular plate-like arrangement. Cells within the arrangement differ in shape from those outside. Outer cells form spine-like extensions of the colony. Cells are also non-motile (no flagella).

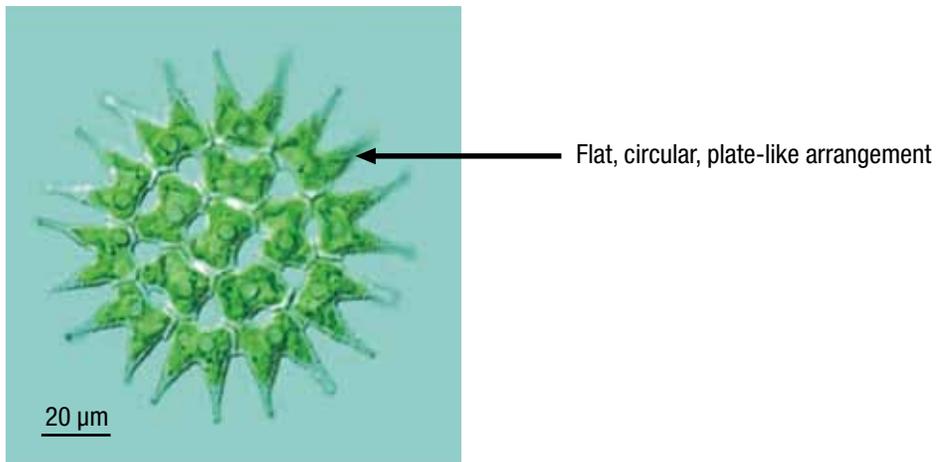

D. Patterson and M. Farmer (micro\*scope 2006)


---

## ► STIGEOCLONIUM

---

A poster-genus for green algae, *Stigeoclonium* is a distinctly bright green alga often attached to rocks in fast-flowing water. Although it forms mats or streamers, it is never as large or as slimy as those of *Spirogyra* (pg30). Colonies will definitely not pass through the open fingers of a cupped hand.

**Under the microscope:** Identification to species are based on branching patterns of the colonies and by the overall shape of the formed tufts. Filaments grow as erect branched tufts or plumes, and the branching can either be alternate or opposite. The branches are thorn-like, short or long, and tapering to fine points.

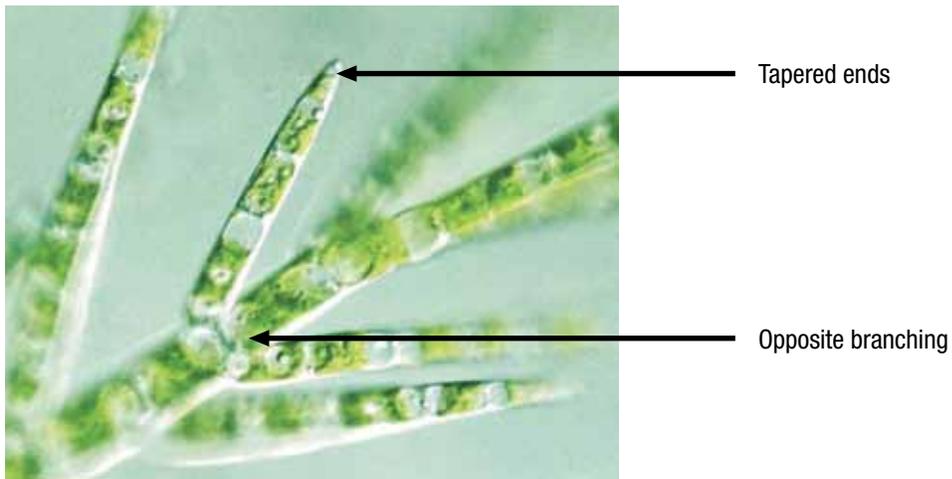


P.V. York (York et al. 2002) (22  $\mu$ m cell width)

## GR Green Algae

### ► OEDOGONIUM

Scooped up from the water, it adheres to one's hands rather than slipping through one's fingers like *Spirogyra* (pg30). *Oedogonium* frequently begins life as small filaments attached underwater to near-shore vegetation. Some species remain attached, and excessive growth can impede water flow in drainage canals. Some detach to form free-floating, fluffy colonies found just beneath the water's surface. Colonies are usually light green in youth, fading to a pale green or tan colour with age. They often occur in such dense associations near shore that evaporating water leaves behind dried masses of 'algal paper'. *Oedogonium* is an unbranched filament that is common in quiet fresh waters. This genus is also very susceptible to copper treatment.

**Under the microscope:** *Oedogonium* cannot be identified well to species if samples are collected before reproductive structures have developed. Also, many species may grow together in close association, making identification to species a challenge. Filaments of cells are not quite cylindrical but are slightly larger at the anterior end.

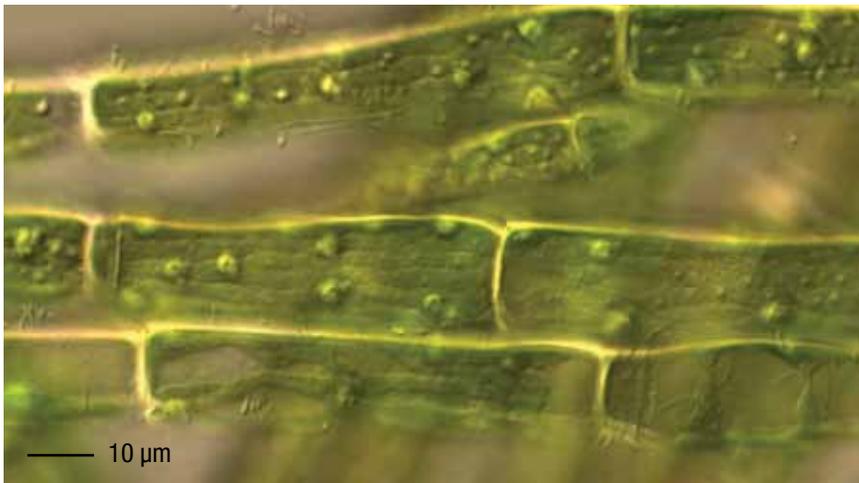

D. Patterson (micro\*scope 2006)

---


## ► CLADOPHORA

---

There are a few freshwater species of this genus. Freshwater *Cladophora* species are commonly found as colonies attached to rocks, plants and almost any other stationary aquatic substrate, but can be detached by wave action. And although they are epiphytic (attaches itself to other substrate or organisms), they produce no mucilage and hence are heavily fouled by epiphytes. They occur in both flowing and standing water, although they are most often associated with streams, rivers and waterfalls.

*Cladophora* growth results in long, filamentous mats. The form of the filaments within the mat can vary between habitats and seasons. Mats that overwinter or constantly wave-washed will have a sturdier form than those found in slower-moving or standing water. *Cladophora* can be fairly abundant during warmer periods of the year, and ranges in color from green to yellow-brown to orange. Septic odour is present when they are abundant, which can be attributed to the decomposition of the algal mat.

**Under the microscope:** Filaments are branching, not tuft, and composed of multinucleate cells. Cells walls are thick, and chloroplasts can be reticulate (net-like) or divided into numerous small discs.

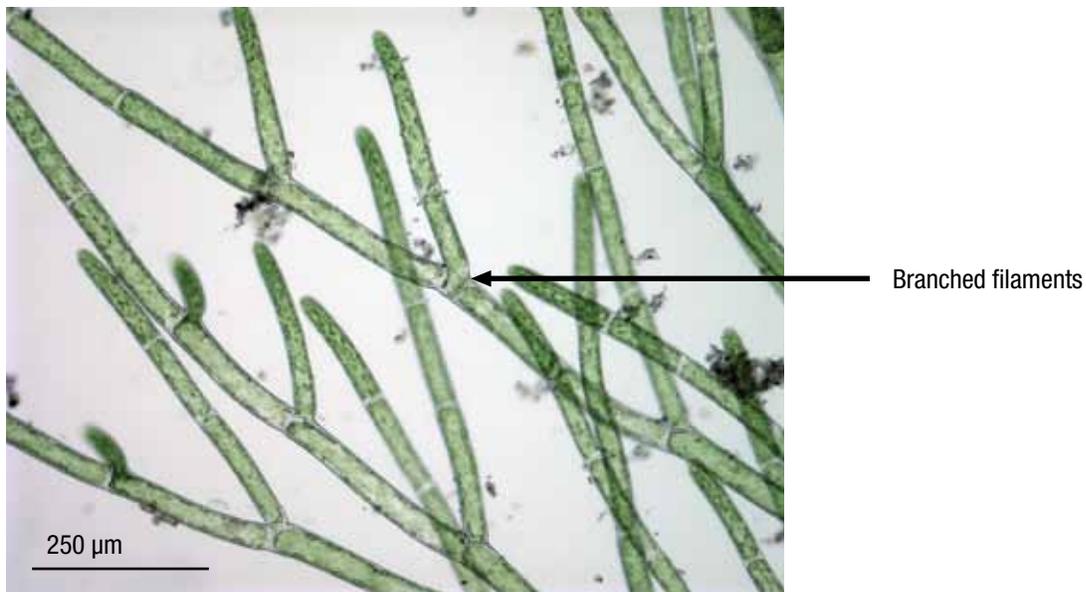

W. Bourland (micro\*scope 2006)

## GR Green Algae


### ► SPIROGYRA

*Spirogyra* or 'water silk' is a common green alga found in shallow warm water where it can form extensive floating mats. It is the most frequently encountered member of the 'frog spit' Order Zygnematales. One can detect the silky quality of the filaments by trying to lift the alga from the water. This genus is very common in both standing and flowing waters. Species can be either attached or free-floating.

The floating varieties can form large mats capable of blocking stream channels and are visible as cloudy filaments just below the water's surface. Colour is usually a vibrant green while the alga is in its spring growth stage. Once it reaches a reproductive stage however, the colour can shift to a yellowish-brown. The filament mats feel very slimy to the touch, noticeably so when swimming or wading through them. Species of *Spirogyra* can be found from mid-summer to fall. Abundant growth in spring can result in clogged filters, seepage areas, drainage canals, and coves. Dense mats can deplete oxygen, interfere with fishing, or provide breeding sites for mosquitoes.

**Under the microscope:** The reasoning for the name *Spirogyra* is evident. While the filaments themselves do not spiral, the chloroplasts within the filament cells certainly do. The chloroplast is a bright grass-green helical peripheral ribbon. Identification to species is not possible without mature specimens. Try not to confuse the name *Spirogyra* with *Spirulina* (a cyanobacteria, pg14). They each spiral in their own special way.

**Under the microscope:** Filaments are branching, not tuft, and composed of multinucleate cells. Cells walls are thick, and chloroplasts can be reticulate(net-like) or divided into numerous small discs.

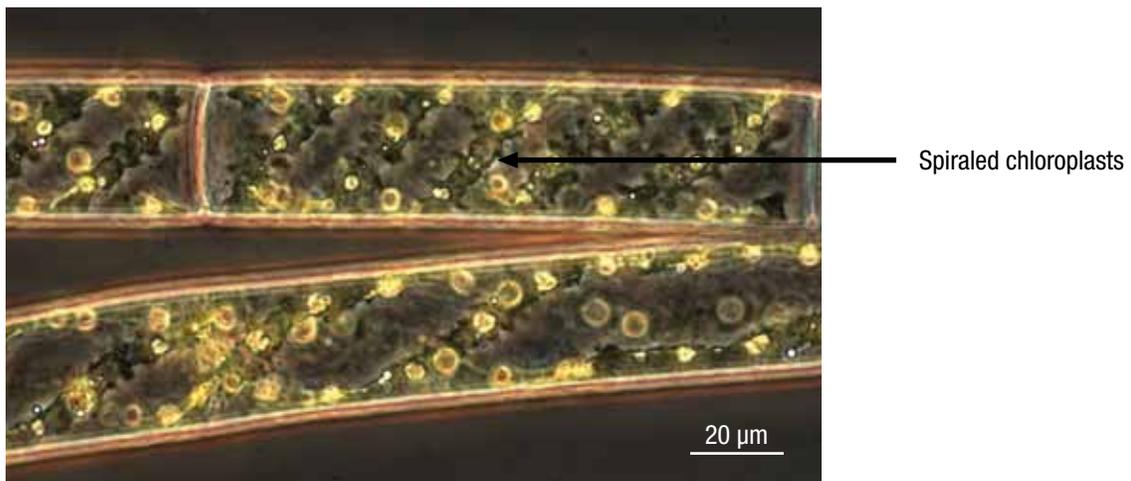

Australian Biological Resources Study (micro\*scope 2006)


---

## ► CLOSTERIUM

---

*Closterium* is generally found free-floating in standing water, and usually mixed in with other species of green algae, especially filamentous forms. It is common in acidic (pH<7) waters, while a few species are found in nutrient-rich hard waters.

**Under the microscope:** Species in this genus are often single-celled (solitary) and may have a variety of shapes but are commonly elongated and crescent-shaped (bowed). They follow a pattern of being slightly widened in the middle and pointed on the ends.

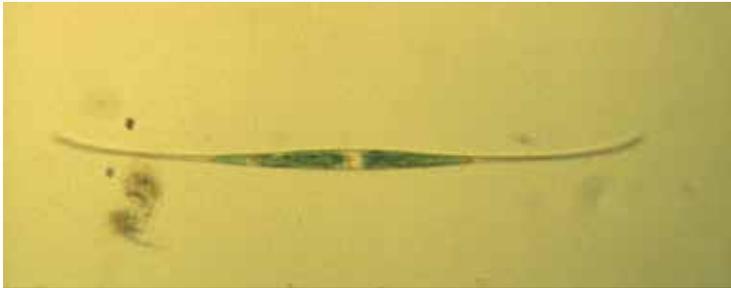

M. Graham and S. Murrell (400X magnification)



D. Patterson and A. Laderman (micro\*scope 2006)


## GR Green Algae

### ► CHARA

*Chara*, or stonewort, is a green alga that is plant-like in appearance and can grow several decimetres tall. Most grow in still, clear freshwater where they form extensive underwater meadows, anchored into the sediment by their rhizoids. They are easily recognized by their whorled branches and musky or garlicky odour.

They are particularly abundant in hard, basic waters (pH>7) where calcium is abundant. Calcareous deposits are found on the surface of the algae, making them feel rough to the touch, which gives rise to its name 'stonewort'. Problems with the species occur when dense growths impede water flow and interfere with recreational activities.

It is often confused with *Nitella* (not referenced in this manual), as both look very similar. Branches arise from each node: with *Nitella*, the branches themselves also branch; in *Chara*, the branches do not branch. Also, *Nitella* is not ill-smelling and found in soft water or acid lakes and hence, is not encrusted with lime like *Chara*.

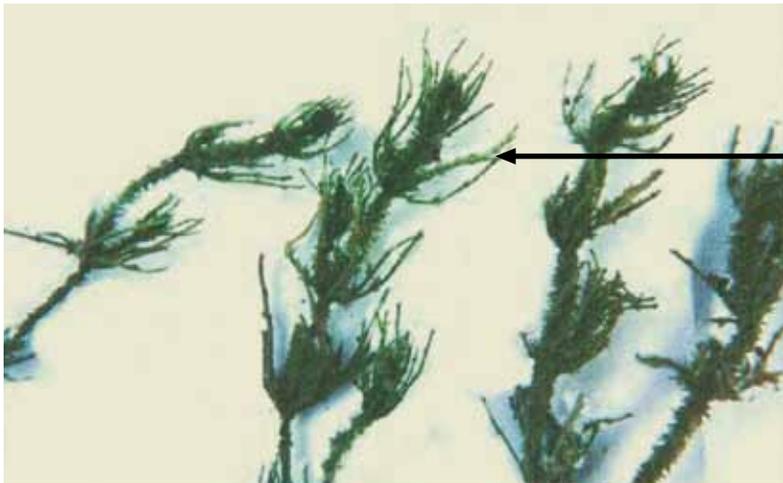

Whorled branches

N.F. Stewart (York et al. 2002) (In situ)


---

**► ASTERIONELLA**

---

Colonies are too small to be visible to the naked eye, but they are planktonic and often very abundant in hard water lakes. They are well known to clog filters because of their rigid walls. An easy way to identify their presence is by smell- if low to moderate numbers exist, a geranium odour will likely be present. However, the odour changes to fishy when high numbers are present.

**Under the microscope:** Frustules (unique cell wall of silica) in girdle view are elongated and rectangular, and form a circular colony (usually 8 cells) in which the cells radiate from a common center like spokes of a wheel.

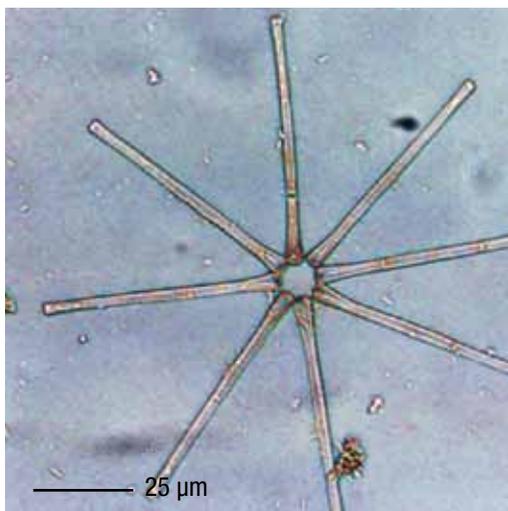

M. Vis (Ohio University 2002)

# DM Diatoms

---

## ► SYNEDRA

---

*Synedra*, in moderate to high abundance, has a distinctive earthy to musty odour. It is also well known to be a typical filter-clogging alga because of its rigid walls. *Synedra* is common in the plankton community and in scums on substrate. Colonies are too small to be visible to the naked eye.

**Under the microscope:** *Synedra* species are large, elongated needle-shaped diatoms that are often bilaterally symmetrical. They are much longer than wide and are usually straight, although some species are curved, in which case they are not bilaterally symmetrical. Frustules (unique cell walls of silica) are either solitary or in radiating colonies attached to substrate at one end.



M. Vis (Ohio University 2002)


---

**► CYCLOTELLA**

---

Besides *Synedra* (pg34), *Cyclotella* is a common bloomer of the diatoms. It is a well-known filter-clogging alga and often found in 'clean' water. That is to say, its abundance often decreases in nutrient-rich water. Colonies are too small to be visible to the naked eye but the species are planktonic.

**Under the microscope:** The valve view of *Cyclotella* will show circular cells with a smooth central region and a grooved or ridged outer edge. The girdle view is narrowly rectangular.

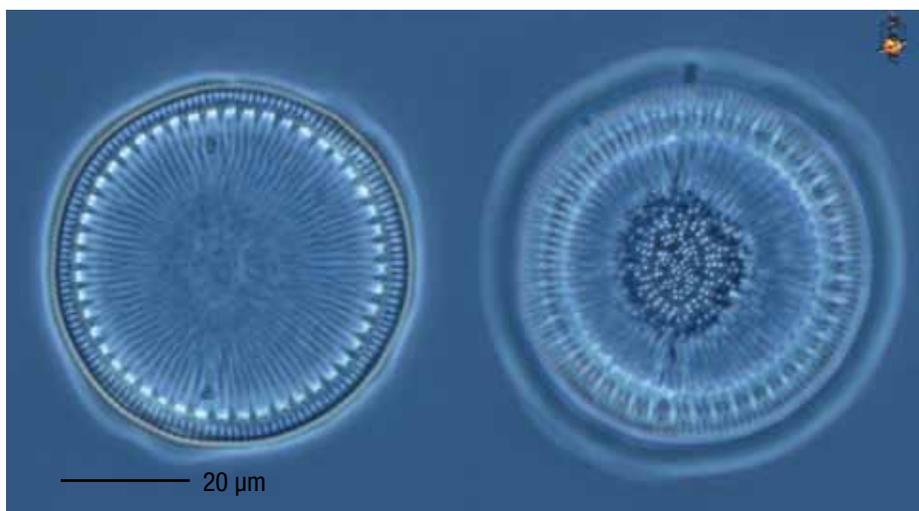


Valve view: M. Bahr and D. Patterson (micro\*scope 2006)

# FL Flagellates

## ► EUGLENA

There are more than 800 species in this genus. They are mostly found in freshwater habitats that are polluted by organic waste or decaying organic matter which can give rise to very dense populations of bloom proportion. *Euglena* blooms can be red.

**Under the microscope:** Euglenids are motile (have flagella), typically spindle-shaped and often uni-cellular. Many have disc-like chloroplasts (pigmented structures within the cell) that are able to supplement photosynthesis by taking up organic compounds. Colourless species do exist however, and are totally dependant on heterotrophic nutrition. Cells lack a distinct cell wall, do not stain (negative test for iodine) and are characterized by their red eye spot.

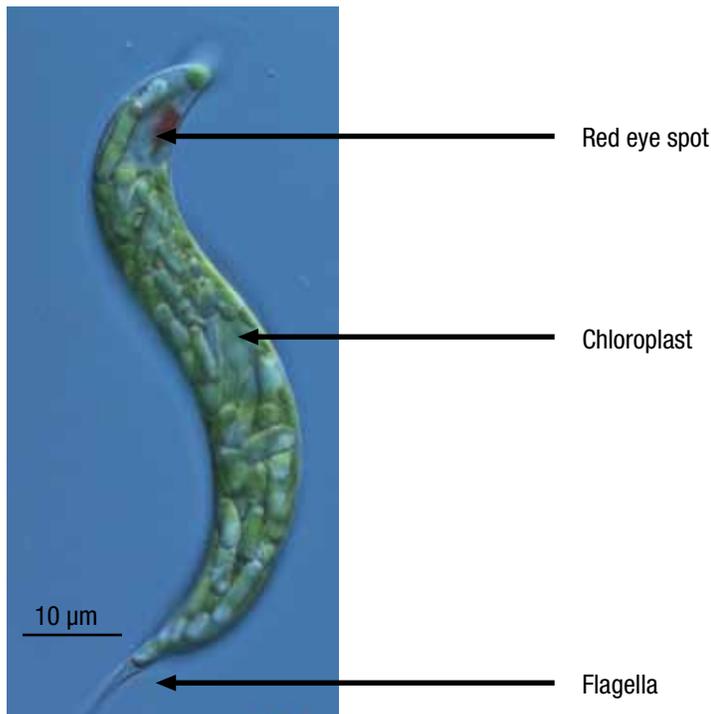

D. Patterson and M. Farmer (micro\*scope 2006)


---

► **SYNURA (GOLDEN ALGAE)**

---

*Synura* species are common in hard water lakes. When present in high abundance, they can cause odour and taste problems. *Synura* exhibits a ripe cucumber/musk melon odour in low to moderate numbers, and a fishy smell in high abundance with other algae. When tasted, it is bitter and leaves a dry sensation on the tongue (tasting algae is not recommended however).

**Under the microscope:** Cells are elongate-ellipsoid and compactly arranged to form a spherical colony that is NOT in a mucilaginous sheath. Species which produce siliceous covering on the surface of the cell may have bristles or scales. Cells are biflagellate, having flagellae of equal length. Chloroplasts are golden brown, a common characteristic of Chrysophytes (golden algae).

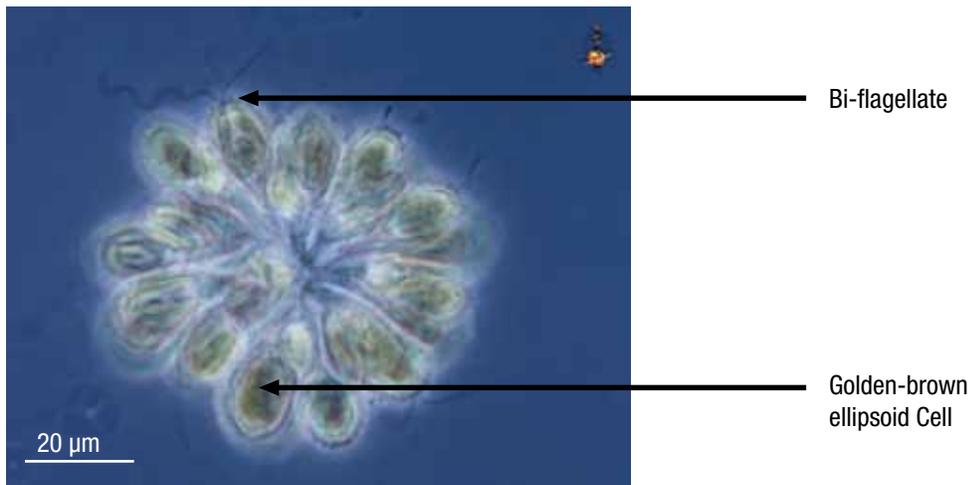


D. Patterson and A. Laderman (micro\*scope 2006)

**FL** Flagellates

► **DINOBYRION (GOLDEN ALGAE)**

In summer, *Dinobryon* species are common bloomers of Golden Algae and are typically found in hard water and eutrophic water bodies. It is an alga of concern because it can cause taste and odour problems, especially in combination with other algae. They are also known to clog filters.

**Under the microscope:** *Dinobryon* species of this genus are characterized by having motile cells enclosed within a colourless envelope. They are colonial and typically have vase-shaped cones (loricas) that are stacked on top of each other which result in a branching chain. Each individual cell has 2 unequal flagellae but the cells are often absent from the lorica. Chloroplasts are golden brown, a common characteristic of Chrysophytes (golden algae).

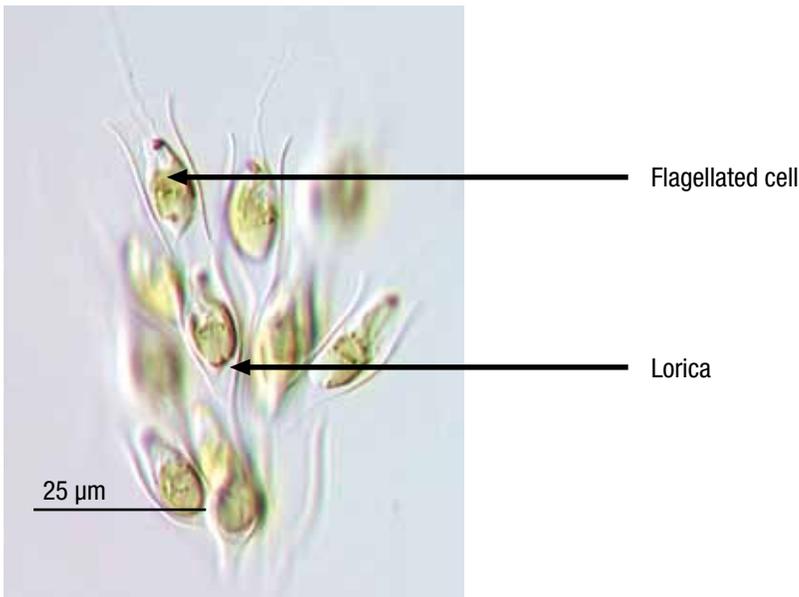

W. Bourland (micro\*scope 2006)

---

**► DUCKWEED**

---

Duckweed is often referred to as algae but it is in fact a tiny colonial, aquatic perennial plant. Individuals are small, disc-like and free-floating. Individuals are typically 2-5 centimeters long. By taking a closer look at the underside of an individual, short hair-like roots will usually be found. Common genera of duckweed include *Lemna*, *Wolffia*, or *Spirodela*.



N. Serediak (in situ)



# Glossary

---

**Akinetes** - a resting spore of cyanobacteria.

---

**Chloroplast** - pigmented structure within the cell that conducts photosynthesis.

---

**Coenobium** - a colony that always has the same number of cells once it is formed.

---

**Epiphytic/epiphyte** - an organism attached to substrate or another organism.

---

**Heterocysts** - a large, thick-walled cell involved in nitrogen fixation in cyanobacteria.

---

**Filamentous algae** - long, multicellular threadlike filaments or masses of algae.

---

**Flagella(e)** - whip-like tail(s) attached to the cell used for locomotion.

---

**Frustules** - unique cell walls of silica.

---

**Lentic** - standing water such as wetlands (ponds, sloughs), dugouts, reservoirs and lakes.

---

**Phytoplankton** - microscopic floating plants, mainly algae that live suspended in bodies of water.

---

**Planktonic algae** - single-celled organisms or aggregations of single-celled algae.

---

**Reticulate** - a net-like pattern; net-like.

---

**Trichome** - a thread of cells.

---

## Photo credits

Microcystis: D. Krogmann and M. Schneegurt. 2006. Purdue University. Cyanosite.  
<<http://www-cyanosite.bio.purdue.edu>>; D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Anabaena: M. Graham and S. Murrell; Australian Biological Resources Study. 2006.  
Micro\*scope. <<http://starcentral.mbl.edu/microscope/portal.php>>.

Aphanizomenon: W. Bourland. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Oscillatoria: B. Anderson and D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Gloeotrichia: P. York in York et al. (2002).

Nostoc: M. Bahr and D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Gloeocapsa: M. Vis. 2002. Ohio University. Algae Homepage.  
<<http://vis-pc.plantbio.ohiou.edu/algaeindex.htm>>.

Chroococcus: D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Chlorella: W. Bourland. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Tetraspo D. Patterson and M. Farmer. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Closterium: M. Graham and S. Murrell; D. Patterson and A. Laderman. 2006.  
Micro\*scope. <<http://starcentral.mbl.edu/microscope/portal.php>>.

Spirogyra: Australian Biological Resources Study. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Oedogonium: D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Hydrodictyon: M. Vis. 2002. Ohio University. Algae Homepage.  
<<http://vis-pc.plantbio.ohiou.edu/algaeindex.htm>>.

Stigeoclonium: P.V. York in York et al. (2002).

Cladophora: W. Bourland. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Chara: N.F. Stewart in York et al. (2002).

Chlamydomonas: D. Patterson and A. Laderman. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Synedra: M. Vis. 2002. Ohio University. Algae Homepage.  
< <http://vis-pc.plantbio.ohiou.edu/algaeindex.htm>>.

Cyclotella: M. Bahr and D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Pediastrum: D. Patterson and M. Farmer. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Scenedesmus: L.A. Zettler and D. Patterson. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Dinobryon: W. Bourland. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Synura: : D. Patterson and A. Laderman. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Euglena: D. Patterson and M. Farmer. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

Duckweed: N. Serediak

Asterionella: M. Vis. 2002. Ohio University. Algae Homepage.  
< <http://vis-pc.plantbio.ohiou.edu/algaeindex.htm>>.

Spirulina: D. Patterson, L.A. Zettler and V. Edgecomb. 2006. Micro\*scope.  
<<http://starcentral.mbl.edu/microscope/portal.php>>.

## References

- Bold, H. and M. Wynne. 1978. Introduction to the algae structure and reproduction, 2ed. Prentice-Hall, INC., Englewood Cliffs, N.J. 662pp.
- Hickman, M. 2000. A brief introduction to algae- Botany 333, Laboratory Manual. University of Alberta, Edm. Canada.
- John, D.M., B.A. Whitton and A.J. Brook. 2002. The Freshwater Algal Flora of the British Isles: An Identification Guide to Freshwater and Terrestrial Algae. The Natural History Museum. Cambridge.
- Palmer, Mervin. 1962. Algae in water supplies: An illustrated manual on the identification, significance, and control of algae in water supplies. U.S. Department of Health, Education and Welfare. Public Health Service. Division of water supply and pollution control. #Pub.-216489.
- Palmer, Mervin. 1977. Algae and water pollution. EPA-600/9-77-036. Environmental Protection Agency, United States of America.
- Prescott, G.W. 1964. How to know the fresh-water algae- an illustrated key for the identifying the more common freshwater algae to genus, with hundreds of species named and pictured and with numerous aids for the study. Dubuque, Iowa, W.C. Brown Co. 293pp.
- Smith, G.M. (ed.), L.R., Blinks, H.C., Bold, K.M., Drew, F., Drouet, J., Feldmann, F.E., Fritsch, H.W., Graham, O.P., Iyengar, D.A., Johansen, H.H., Johnson, B.H., Kethcum, G.F., Papenfuss, G.W., Prescott, E.G., Pringshsheim, H.H Strain, and L.H., Tiffany. 1951. Manual of Phycology: An introduction to the algae and their biology. Chronica Botanica Co. (Pub). Waltham, Mass., USA.
- P.V. York, D.M. John, and L.R. Johnson. 2002. A Photo Catalogue of Images of Algae and Algal Habitats. A CD-ROM accompanying 'The Freshwater Algal Flora of the British Isles' (eds. D.M. John, B.A. Whitton, A.J. Brook). Cambridge University Press.

## Recommended reading

### ***Algal treatment and removal, and management practices***

Alberta Agriculture, Food and Rural Development, Agriculture Education and Training Branch. 2002. Quality Farm Dugouts. Prairie Water News. Alberta, Canada.

Palmer, Mervin. 1962. Algae in water supplies: An illustrated manual on the identification, significance, and control of algae in water supplies. U.S. Department of Health, Education and Welfare. Public Health Service. Division of water supply and pollution control. #Pub.-216489.

### ***Field sampling methods***

John, D.M., B.A. Whitton and A.J. Brook. 2002. The Freshwater Algal Flora of the British Isles: An Identification Guide to Freshwater and Terrestrial Algae. The Natural History Museum. Cambridge.

### ***Preserving and preparing mounts***

Prescott, G.W. 1964. How to know the fresh-water algae- an illustrated key for the identifying the more common freshwater algae to genus, with hundreds of species named and pictured and with numerous aids for the study. Dubuque, Iowa, W.C. Brown Co. 293pp.

## On-line resources and identification keys

Algae Vision. 2006. Natural History Museum of London.  
<[Http://www.nhm.ac.uk/index.html](http://www.nhm.ac.uk/index.html)>.

Key to common freshwater algae. Connecticut College. Botany Department. <[http://aspen.conncoll.edu/SiverWeb/LucidKeys/Carolina32\\_Key.html](http://aspen.conncoll.edu/SiverWeb/LucidKeys/Carolina32_Key.html)>.

Micro\*scope. 2006. <<http://starcentral.mbl.edu/microscope/portal.php>>.



